


PERIODONTAL

HOE SCALER

Used to dislodge heavy supramarginal calculus.

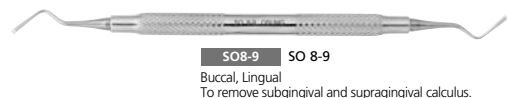

SCALER

The handle is coated with pressure sensitive silicone and provides unique tactile sensitivity, control and color identification.




U15 has durable and strong blade and it is designed to remove of heavy calculus for buccal, lingual and proximal.



FILE

Used to crush large calculus deposits so that the deposit can be more readily removed by a curette. Can also be used to smooth the margins of amalgam restorations.

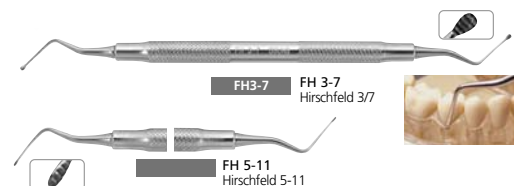

For interproximal To crush and remove heavy deposits from subgingival and supragingival interproximal areas.



STANDARD GRACEY CURETTE


American Type

Curettes are used supragingivally and subgingivally on enamel and root surfaces. There are two major sub-types of curettes: universal curettes and area-specific curettes. Universal curettes can be used to remove medium-sized and small calculus deposits. Area-specific curettes are designed to remove light calculus and endotoxins from root surface. Gracey Curette is an area-specific curette. It has Offset blade that provides a perfect working angulation for the teeth surface.



UNIVERSAL CURETTE

American Type

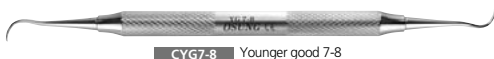

SUB-ZERO CURETTE

