


## RIGID GRACEY CURETTE American Type

The Rigid Curette has a stronger and less flexible shank than the standard one. So it can be used for heavy calculus removal without using additional scalers.



## RIGID UNIVERSAL CURETTE American Type

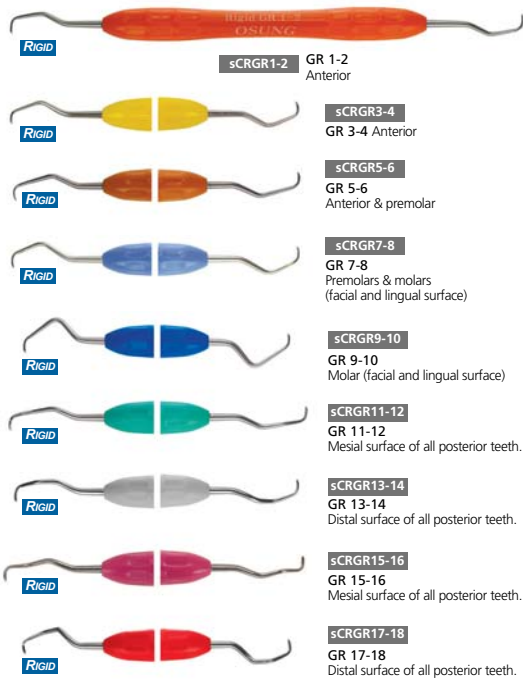


**Gracey Curette Chart**

- GR1-2, GR3-4, GR5-6
- GR7-8, GR9-10
- GR11-12, GR15-16
- GR13-14, GR17-18

**Rigid?**  
Rigid has a strong terminal shank.


## MINI GRACEY CURETTE American Type

The Mini Curette is very useful for accessing into deep periodontal pockets and root surface of 5mm or more.

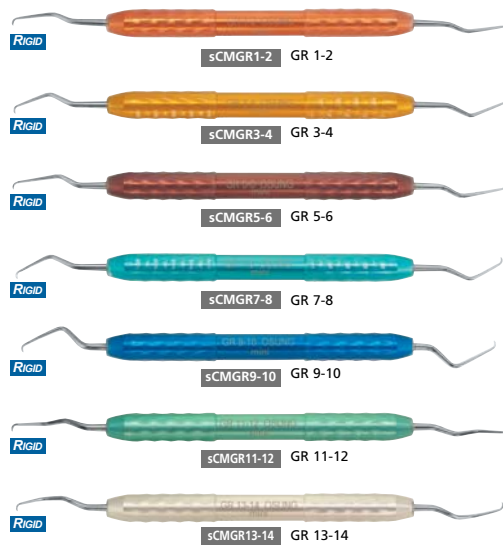

### MINI CURETTE ?


Terminal shank elongated by 3mm for accessing into deep periodontal pockets and root surfaces of 5mm or more. The blade length is reduced in half from the standard Gracey curette, for a better adaptation in narrow pockets and furcations. Blade is a little thinner than the standard to enable easy gingival insertion and reduce tissue distention.

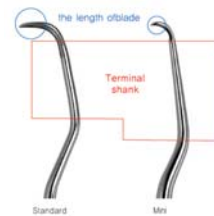

## RESTORATIVE

### COMPOSITE INSTRUMENTS

Non-stick-instruments for a fast and efficient placement and shaping of composite, and an easy clean up. The tips are highly polished and titanium-coated.

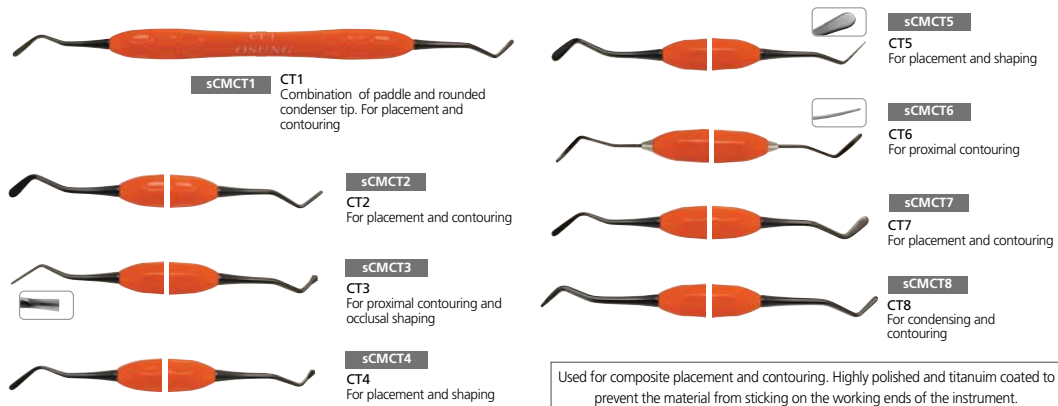

Used for composite placement and contouring. Highly polished and titanium coated to prevent the material from sticking on the working ends of the instrument.